


Matrix LT Air

The Matrix LT Air Mattress system provides alternating pressure therapy with micro low air loss to aid patients who are at risk of developing pressure ulcers or have existing stage 1, 2, 3, or 4 ulcers. Caregiver friendly, the Matrix LT Air is an affordable, therapeutic surface solution for any clinical care setting.

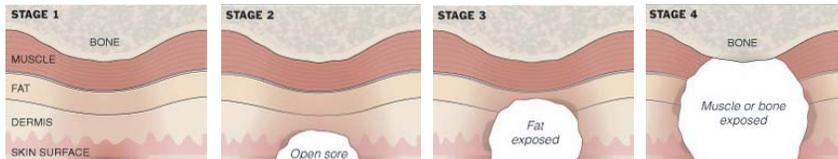

Various top covers are designed and available to satisfy your individual requirement in terms of therapy, prevention and patient comfort.

Master Control Unit

- Dimension: (L)13.5"x (W)7.3"x (H)8.3"
- Weight: 7.1 lbs
- Rating Voltage: 110-120 volts
- Max Current: 0.2A max.
- Rating Frequency: 60Hz
- Cycle Time: 5,10,15,20 mins
- Static Function
- 3 in 1 alternation
- Min/Max Pressure: 12~47 mmHg +/- 6 mmHg
- Air flow rate: 10 liters per minute
- Protection Type: Class I, Type BF
- Low Pressure Alarm: Visual / Audio
- Power Failure Alarm: Visual / Audio
- Maintenance service LED indicator

Mattress Replacement

- Model No: Matrix LT Air
- Dimension: (L)80"x (W)36"x (H)8"
- Weight: 25.3 lbs
- Weight Capacity: 500 lbs
- Cell Material: Blue Nylon with PU Coating
- Cell in cell design to stay inflated during power outage
- Cell Number: 18 cells
- Cover Material: 4 Way Stretch Urethane
- Base Material: Black Polyester with PVC Coating
- CPR Quick Release System
- Cable Management System on bottom cover

Custom Medical Solutions

(866) 350-5640

www.custommedicalsolutions.com

CUSTOM MEDICAL
SOLUTIONS

"Delivering BETTER SOLUTIONS today
for a BRIGHTER TOMORROW!"



Matrix LT Air


HAPPYHEEL™ Heel Relief
Patent Pending



Matrix LT Air


provides you, the caregiver, with peace of mind
and your patient, with world class care & comfort

Perfect for patients and caregivers

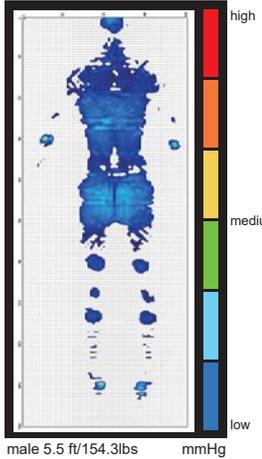

PRESSURE MAPPING
Matrix LT Air provides constant low pressure by redistributing the patient's profile over a greater area than a standard mattress.

Nylon covered bottom cover for friction and noise reduction when raising bed's head pan

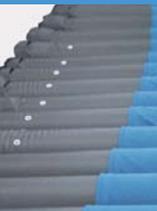

Effective micro low air loss for moisture management



Slimline pump design is unobstructive

8 adjustable comfort settings provide different levels of comfort depending on the patients therapy requirement.

HAPPYHEEL™ provides further therapy and comfort by decreasing pressure in patient's heel area.



- Interchangeable air cell design provides easy cleaning and minimizes the proliferation of bacteria
- The Matrix LT Air employs 8" air bladders with cell in cell design (4.5" of therapy, 3.5" of static air) to provide therapeutic alternating pressure and also allow the mattress to stay inflated in a power outage.
- The Matrix LT Air replacement mattress can accommodate patient weight up to 500 lbs.

HAPPYHEEL™ provides extra softness in the last 5 air cells to provide protection for patient's fragile heels

• Modular design ensures ease of use and cleaning


Matrix LT Air is your solution for prevention and treatment of pressure ulcers



Cable management system keeps electrical cables off floor for caregiver and patient safety



CPR quick release system deflates the system within 20 seconds

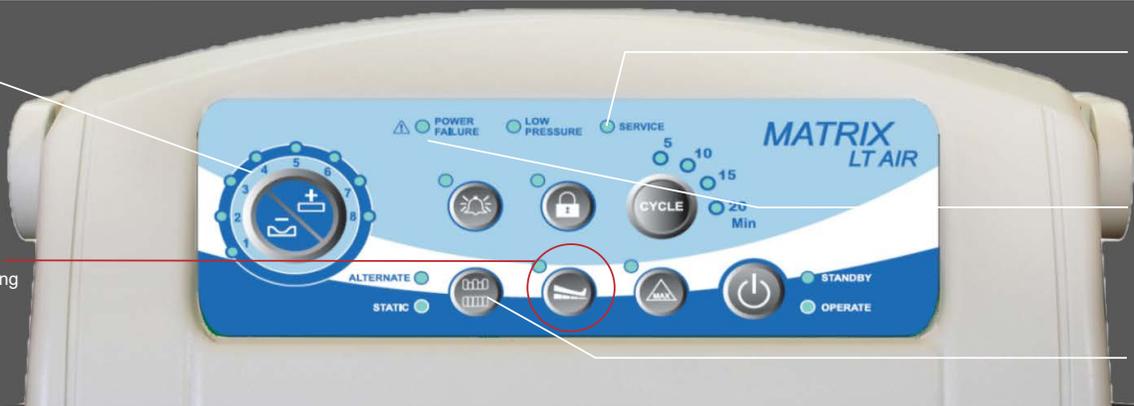

Service LED and alarm indicator is a friendly reminder for advising caregiver to send back control unit to certified technician for annual inspection.

The Power Failure alarm warns in case of power outage.

3-in-1 alternating pressure covers 66% of patient's body (vs. 50% on 2 in 1 systems) to provide a more comfortable yet still effective therapy



Spring loaded hooks snap to fit your bed's footboard